

## Guarantee Certificate

Product Model:

Purchase date:

Purchase place:

Sales Dept:

### One Year Repair Guarantee

1. In the unlikely even that your product develops a fault within the first one years of purchase :  
(Please return your product along with your receipt and the original guarantee certificate and we will repair it free of charge provided you have operated it according to the operating

2. What is not covered?

(1) Not following the operating instructions

(2) Incorrect or unauthorized repair

(3) Accidental damages during your own transportation , or you drop it on the floor or you deliberately damage or neglect the product after purchase .


(4) Fire , earthquake , weather such as flood, lightning , high wind.

3. Please keep your guarantee certificate safe, as no duplicated certificate will be issued.

## Ultrasonic Fog Moisturizing Beauty Instrument

## Ultrasonic Fog Moisturizing Beauty Instrument

### User's Manual

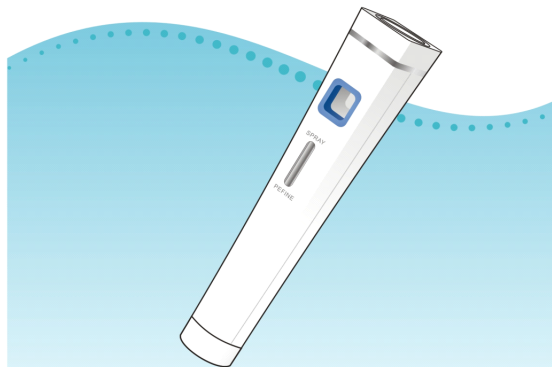

- Thank you for selecting our product.
- Please carefully read the User's Manual for normal operation before use the product. Please keep it for future read.
- This beauty instrument is for healthy skin treatment with the purpose of beautifying skin. It is not allowed to use or sell it for therapy purpose.

## Contents

I. Applicable Skin Defects	1
II. Innovative Technology	1
III. Safety Cautions	2
IV. Introduction to Parts	4
V. Usage	6
VI. After Use	11

## Specification

Name: Ultrasonic Fog Moisturizing Beauty Instrument

Power Supply: Three 1.5V7# batteries

Dimension: (Body) 17 (H) X 3.8 (W) X 4.3 (T) cm

Weight: (Body) 178g

## I. Applicable Skin Defects

Temperature change caused by air conditioning, environmental pollution and season change, pressure, ageing and slow metabolism may cause water loss in skin, and result in rough and dull skin, dry lines and diversified skin defects.

- ★ **Mature Skin:** sufficient water in skin may moisturize the mature and ageing cells which have slow metabolism, recover skin to smooth state, fight ageing and prevent wrinkles.
- ★ **Skin with Acne:** regularize hormone secretion, replenish water, resist external irritation and prevent acne.
- ★ **Sensitive Skin:** maintain sufficient water in skin cells to prevent dull and lifeless tone.
- ★ **Rough Skin:** constantly moisturize skin to keep subcutaneous fat at "semi liquid state", shrink rough pores, and make skin moist, soft and firm.

## II. Innovative Technology

This instrument adopts nanometer technology and ultrasonic technology and can spray nourishing skin toner into nanometer mist onto skin, massage the cells with the ultrasonic massage top which vibrates for 3 million times/second, and can hydrate and nourish skin, assist to firm fine lines, and make skin smooth and tender.

### III. Safety Cautions

- This instrument is special for beauty care for healthy skin, please do not use it for other purposes.
- Please follow the instructions on use frequency and operation, please shorten the duration when use for the first time or use it for sensitive skin.
- Please stop use if you feel uncomfortable when use.

#### Forbidden users and parts:

- It is not allowed to use for forehead, temple, eyeball, throat and bone joints; Otherwise, it may possible cause skin or physically uncomfortable;
- Skin which is hot or inflammatory after exposing under sunlight;
- Dermatitis caused by cosmetics;
- Part with acne;
- Spots which feel itching or hot;
- Spots caused by physical irritation;
- Users with angioedema caused by steroid, hormone and liver dysfunction;
- Face-lifting part;
- Do not use the instrument for the same part on the same day after depilating and shaving;
- Other parts which are not specified in the User's Manual.



#### Cautions on Operation and Storage:

- Do not use it when the power cord, plug and socket are damaged;
- Do not store it on washing table which is moist;
- Do not operate it with a moist hand;
- Do not store it under power on state;
- Keep it away from children, Do not use it for minors.

- Please perform dermal sensitivity test before use for the first time

#### Simple method for dermal sensitivity test:

1. Clean the wrist;
2. Extend Special Essence onto wrist for 48 hours, dry it naturally (do not wet skin when test, do not rinse off the essence), immediately stop the test when skin feels uncomfortable;
3. You may use the instrument if skin doesn't feel uncomfortable within 48 hours;
4. Please immediately stop use and see a doctor when allergy occurs



- 2 -



### Important safety notice



Do not store it in the environment with high temperature, keep far from fire source!



Keep it away from children!



It is not allowed to disassemble the instrument by yourself



Please remove the batteries from the instrument when it won't be used for a long period




Do not use it in water, avoid dampness!




It is not allowed to drop and knock the instrument!

- 3 -

### IV. Introduction to Parts

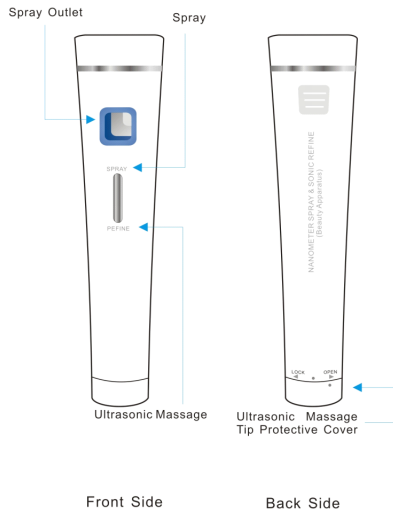

- 4 -

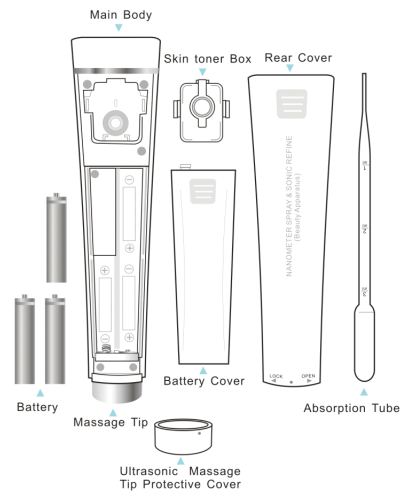

- 5 -

**Ultrasonic Fog Moisturizing Beauty Instrument has two functions**


- 1 Replenish nutrition and water with nanometer mist;
- 2 Hydrate and nourish skin, lift and firm skin with ultrasonic massage.

**V. Usage**


**Preparation:**

- Install the batteries and skin toner --- please use chargeable batteries for environmental protection

- 1) Remove the protective cover of the Massage Tip

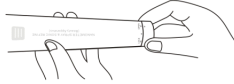


- 2) Open the rear cover

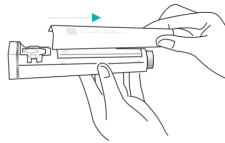


- 6 -


- 3) Remove the battery cover

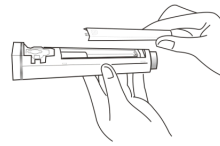

- 4) Install batteries

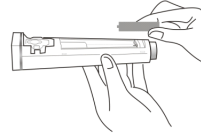

- 5) Fill into the skin toner into the small box, install the rear cover.

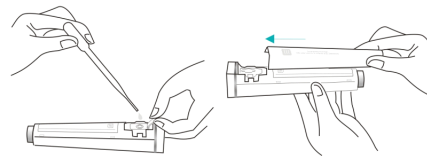

- 7 -


**● Step 1: Replenish nutrition and water with nanometer mist**

Press the Spray key; start nanometer mist, skin toner sprays onto the face as nanometer fog to moisturize skin



**● Step 2 Hydrate and nourish skin, lift and firm skin with ultrasonic massage**

1. Remove the protective cover of the Massage Tip and put it onto the top cover;

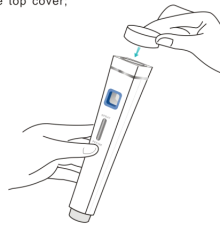


- 8 -

2. Extend gel onto face and massage tip;




3. Press the REFINE key to start ultrasonic massage and repair skin.

1. Operate the "Ultrasonic Fog Moisturizing Beauty Instrument" around the chin. Attach, press and push upwards along the curve of cheek, and slowly move to skin behind the ear.

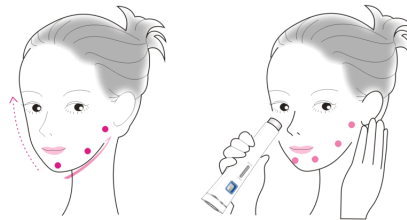

- 9 -

2. Attach, press and push downwards along the neck, and slowly move to the shoulder to promote lymph fluid flowing from the skin below eartowards to the collarbone.

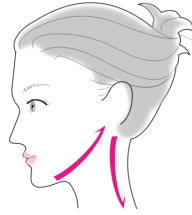

3. Attach the "Ultrasonic Fog Moisturizing Beauty Instrument" onto skin and push and pull upwards according to the drawing below.



4. Repeat for 5-10 times for every above step.

-10-


## VI. After Use

After Use: clean the instrument, remove the batteries, and store it in the protective sleeve.


★ Clean the skin toner box:

- After using up the skin toner, please clean it according to the following steps:

1. Remove the protective cover of the massage tip;

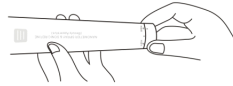


2. Remove the rear cover;

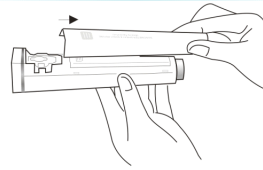


-11-


3. Remove the skin toner box;

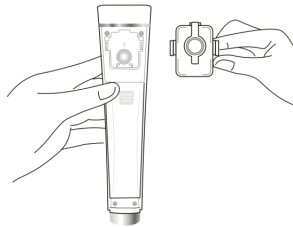

4. Clean it and then install into the instrument,

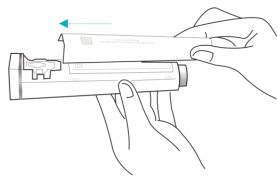

-12-

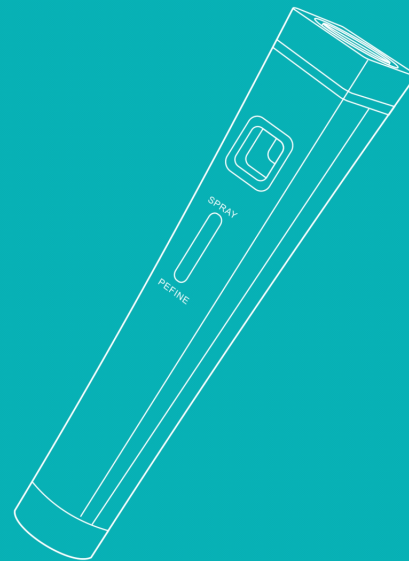

Ultrasonic Fog Moisturizing Beauty Instrument