

# ULTRASONIC-IONIC REPAIRING INSTRUMENT

## 【 Quality Warranty Card 】

Product Model

Purchased On

Purchased At


Distributors

- The warranty period of this product is 1 year

1. Please use the product according to the Operation Manual, faults within the warranty period enjoys free repair.  
(Please provide the product and the Warranty Card for free repair service)
2. We provide repair service with certain charge if one of the followings occur even within the warranty period.
  - (1) Fault and damage caused by incorrect use;
  - (2) Fault and damage caused by unauthorized repair or reconstruction;
  - (3) Fault and damage caused by transportation and falling down after purchasing the product;
  - (4) Fault and damage caused by fire, earthquake, flood, lightning strike and other force majeure.
3. Please carefully keep the Warranty Card since it shall not be reissued.

## ULTRASONIC-IONIC REPAIRING INSTRUMENT

### Operation Manual

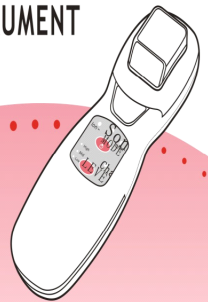

- Thank you for selecting our product.  
Please carefully read the Operation Manual before use for correct operation.  
Please take good care of it after reading.
- This beauty instrument is based on the treatment for healthy skin with the purpose of beauty care, it is not allowed to use or sell the product for the purpose of therapy.

## Contents

### Contents

Safety Cautions .....	1
Specification .....	
Introduction to Parts .....	2
Deep cleansing .....	
Induct-in .....	
Ultrasonic and Ionic work together .....	3
Warranty Card .....	4

### Specification

Power Supply: 3 pcs AAA normal or recharge batteries  
Spec. : (Body) H16× W4.A× T 3.5cm  
Weight: (Body) 77.4g  
Frequency of Ultrasonic: 3±5% M Hz

## Safety Cautions

- The product is the beauty care instrument for healthy skin, please do not use it for other purposes except beauty care.
- Please follow the use frequency and operation manual, please shorten the use time for use for the first time and allergic skin.
- Please stop using if the user feels uncomfortable.

### Users and parts which are not allowed to use the instrument

• Eyeballs;  
• Users with fracture repair surgery, face-lifting, metal, silica gel and plastic implanted in face;  
• Hot and inflammatory parts caused by intensive sunlight;  
• Dermatitis caused by cosmetic;  
• Acne and inflammatory parts;  
• Itching or hot spots;  
• Spots caused by physical irritation;  
• Users with angioedema caused by long time use of steroid, hormone and liver dysfunction;  
• Plastic parts;  
• Users with metal in body after therapy;  
• Users with heart disease (particularly for users with heart pacemaker);  
• Other parts which are not designated in the Operation Manual.

### Cautions on Operation and Storage

Do not put it on the washstand and other places which may easily suffer from water;  
Do not operate with wet hand;  
Cut off power supply before storage;  
Keep it out of reach of children, do not use it for young persons.  
★Please perform skin allergy test before use it for the first time

### Simple Method for Skin Allergy Test:

1. Cleanse the wrist;
2. Extend the Special Essence onto skin for 48 hours, keep it dry naturally (skin shall not contact water, do not rinse off the Essence), immediately stop test when skin feels uncomfortable;
3. Use it if skin doesn't feel uncomfortable within 48 hours;
4. Immediately stop using and see a doctor when skin becomes allergic.



## Instrument Body

Face Electrode

Cosmetic Cotton Cover


Cosmetic Cotton Open button

Operation Panel

Hand Electrode

## Parts Introduction

## Introduction to Names on Operation Panel

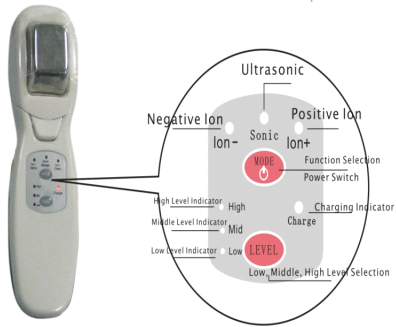

— 2 —

## How to hold hand electrode?

For normal way, you should  
Tightly hold the handle of the  
Beauty Instrument, and which shall only  
be operated by yourself.



— 3 —

Insert the electrode wire to the  
special hole  
and stick the pad on your skin.  
Thus you can care your self  
without holding the hand  
electrode.  
You can also do caring for your  
friend with pad on your friend's  
skin.



## How to power machine and recharge machine.



Machine can be  
powered by 3 pcs  
normal AAA  
battery



Machine can be  
powered  
by 3 pcs recharge AAA  
battery  
with next two  
recharge solution:

1. Recharge with your  
computer USB source  
by connecting  
machine with  
computer with USB  
line.



When recharge,  
the charge light on.

2. Recharge with the adaptor by insert  
adaptor plug into wall socket.



— 4 —

## Deep Cleansing

1. Makeup Removing



## Remove ornaments

Assure to remove earring and  
ear stud, rings, necklace and  
watch.

## Operation instruction

I. Wet the cosmetic cotton with cleaning toner or distilled water  
(make the face and cosmetic cotton damp)

II. Fix the cosmetic cotton on the massage tip

III. Start the instrument, select Ion+

IV. Select High/Low Level:


Please select LOW when use for the first time, and then select  
Middle or High after skin has adapted to it

## V. Treatment Method

Attach the electrode Pad onto skin, or tightly hold the handle,  
attach the massage tip onto face, move it once every 10-30s, the  
moving direction is shown in the right diagram.

## VI. After Treatment

Open the cosmetic cotton cover, remove the cosmetic cotton  
and clean the massage tip

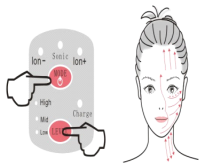


— 5 —

## Induct-in

For the desired part on the face, conduct with a proper amount of essence (for example: essence with spots  
removing or wrinkles removing functions) to increase the concentration and improve the absorption of  
essence.

## 1. Wet the massage tip with essence

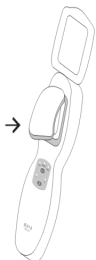

a. Clamp the cosmetic cotton



b. Drop essence onto  
cosmetic cotton -  
thoroughly wet it

each time use  
3-5 g essence

— 6 —



c. Finish

## 2. Start the instrument, select Ion-

## 3. Select High/Low Level:


Please select LOW when use for the first time,  
and then select Middle or High after skin has adapted it

## 4. Treatment Method

a. Extend essence onto desired part, attach the massage  
tip on the part, and keep it onto skin for 20s or longer  
b. Move it once every 5s for other parts (after hearing the  
prompt sound)

## 5. After Treatment

After use, Open the cosmetic cotton cover, remove the  
cosmetic cotton and clean the massage tip

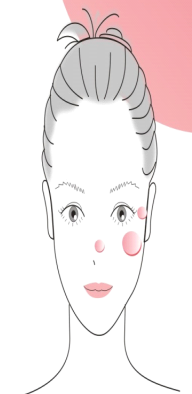

Nursing picture

— 7 —

## Ultrasonic and Ionic work together

Massage the entire face with ultrasonic and Ionic together. Ionic can induct-in nutrition and 3 million times/second ultrasonic can stimulate and adjust the cell membrane, improve the permeability of epidermis, accelerate metabolism and promote skin to absorb essences.



each time use 3-5 g essence




The massage tip shall tightly attach onto skin

2. Start the instrument, select Sonic and Ion-(two indicators light on)

3. Select High/Low Level:

Please select LOW when use for the first time, and then select Middle or High after skin has adapted to it

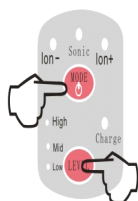

— 8 —

## 4. Treatment Method

“Slim Face”:



operate from mouth to eyes

“Lift Corner of Mouth”:



operate upwards surrounding the corner of mouth

Remove Dark Circle, Pouch, Face Lymph Detoxification

- Operate from inner canthus to temple
- Operate from edge of face to skin behind ear
- From skin behind ear to collarbone part
- From collarbone part to armpit
- It is better to repeat for 20 times

“Firm Double Chin”:



operate from chin to ear, from chin to corner of mouth upwards

T-zone:



slide on nose wing and nose bridge upwards and downwards for 1-2 minutes



5. After Treatment

— 9 —