

CM-5

Photon ultrasonic anti-wrinkle machine

Operating Instructions

Manual



- Thank you very much for choosing this product.
- Please read through carefully and understand these instructions before use.
- Please keep it in a safe place for future reference.
- This device is intended for home use and cosmetic purpose only, using this device for treating any medical condition or illness is not supported by the manufacturer.



Photon ultrasonic anti-wrinkle machine

Photon ultrasonic anti-wrinkle machine



The principal of photon ultrasonic anti-wrinkle skin care machine

- ◆ The ultrasound generator of machine care machine can have 3 million times of micro-amplitude vibration per second. It is most efficient and security ultrasound for face beauty nowday. As it can reach very small distance under the skin, so it is more suitable for the skin care around eyes. The high ultrasonic vibration can give the skin cells moderate stimulation and effectively help activate cells for metabolism.
- ◆ The irradiation with light with different wavelength has different effect to skin: The special pure blue light has the function of calming and against skin rashes; suitable for the skin which is sensitive, oily and with skin rashes. It can adjust skin sebum secretion, effectively get rid of acne and whelk to reach the aim of beautification. Besides, the blue light can astringe and tighten the loosen skin.

The special pure red light is suitable for any type of skin; it can reach derma to stimulate fibre cells. Heating can improve blood circulation and collagen albumen hyperplasia to reduce wrinkles and make skin more elastic and smooth. The mild wavelength is suitable for surrounding eyes and forehead.

The special pure green light between red and blue light is suitable for the mixed skin, it can improve the function of alternating use of cells and oxygen and activate the micro-circulation function of body surface. Meanwhile, calming the skin can make both body and skin relax, so to reduce the pain of skin.
- ◆ The colored light accelerate absorbing the effective element of oily nurture.
- ◆ The over lapping function of tricolored red, green and blue, Galvanic with 3 M ultrasonic waves makes skin care more remarkable.
- ◆ The design of the ultrasonic wave generator and photon care generator keeps close to face and around eyes, thus it is easy to operate.

Contents

◆ The principal of photon ultrasonic anti-wrinkle skin care machine	1
◆ Safety Instructions	2
◆ Parts introduction	4
◆ How to care your face	6
1. daily caring	7
2. Remove Wrinkle Nursing	8
3. Dilute the spot, even skin color, whitening	10
4. Acne caring	10
5. Sensitive skin	11
◆ Guarantee card	12

Specification

Product name: Photon ultrasonic anti-wrinkle machine
Model: ZHF-CM-5
Adaptor input power : AC 100-240V
Adaptor output power : 15V, 400 mA
Dimensions: 170X46X50(mm)
weight : 132g
Ultrasonic frequency : 3MHz±5%

Safety Cautions

- *The product is the beauty care instrument for healthy skin, please do not use it for other purposes except beauty care.
- *Please follow the use frequency and operation manual, please shorten the use time for use for the first time and allergic skin.
- *Please stop using if the user feels uncomfortable.

Users and parts which are not allowed to use the instrument

- ◇ Pregnant;
- ◇ Throat, Eyeballs;
- ◇ Users with fracture repair surgery, metal, silica gel and plastic implanted in face;
- ◇ Hot and inflammatory parts caused by intensive sunlight;
- ◇ Dermatitis caused by cosmetic;
- ◇ Inflammatory parts;
- ◇ Itching or hot spots;
- ◇ Spots caused by physical irritation;
- ◇ Users with angioedema caused by long time use of steroid, hormone and liver dysfunction;
- ◇ Users with metal in body after therapy;
- ◇ Users with heart disease (particularly for users with heart pacemaker);
- ◇ Other parts which are not designated in the Operation Manual.

Safety Instructions

Note: Please read following information carefully before using your machine. This product is well-designed precision instrument. In order to have you device functioning perfectly, please read the below information thoroughly and follow all the instructions. This will help you to ensure the warranty validity and to extend



Simple Method for Skin Allergy Test:

1. Cleanse the wrist;
2. Extend the Special Essence onto skin for 48 hours, keep it dry naturally (skin shall not contact water, do not rinse off the Essence), immediately stop test when skin feels uncomfortable;
3. Use it if skin doesn't feel uncomfortable within 48 hours;
4. Immediately stop using and consult a doctor when skin becomes allergic.

2

Operations and Maintenance Precautions

1. Store your product away from the humidity and do not use in under water. Otherwise the electronic components of your product be corroded or have short circuit.
2. Never attempt to disassemble the product yourself, All the product components can only be maintained and repaired by professional technicians.
3. Don't use any adaptor other than the one included in the package. Using incompatible adaptor may result in malfunction and will void the warranty validity.
4. Keep this product away from children, Children shall not play with this product as the disassembled components may be swallowed by the children accidentally.



Important safety notice



Disassemble the product by yourself is not permitted.



Keep this product away from children




Please keep away from fire.



Store your product away from the humidity and do not use it under water.



Please don't drop and collide it.

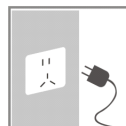

Please take away electricity after usage.

3

Parts introduction

- * Ultrasound with 3 million vibration every second
- * Galvanic
- * High quality stainless steel probe



Special pure red light: 12 pcs LED.
Special pure Blue light: 12 pcs LED.
Special pure Green light: 12 pcs LED.
High energy, fully focused

Galvanic Electrode

4

- Massage intensity indicator
- Ion+ or Ionic-working indicator
- Ultrasound working indicator
- Power indicator
- Massage intensity select
- Ion+ or Ionic- select
- Ultrasound select
- Photon select
- Power on/off select

HIGH
LOW
ION+
ION-
SONIC
POWER
LEVEL
ION+/ION-
SONIC
PHOTON
ON/OFF



• Should touch this electrode when you use Galvanic function

5


Get ready for the use: Operation Procedure and testing

- Place the appliance on a flat surface and put several drops of water on the ultrasonic waves generator.
- Turn NO/OFF on and then press Sonic button. Select one of the levels -Low or high.
- If the water on the ultrasonic waves generator vibrates, the ultrasonic waves are produced.

Suggestion for skincare:

Contact the generator side. (The generator side must tip lightly on your skin)

Put essence on your skin and the probe. Add some essence if the generator side is not lubricating.
(The liquid dosage amount is 4-5g each time.)

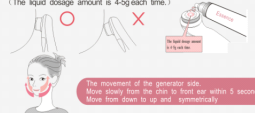

The movement of the generator side:
Move slowly from the chin to front ear within 5 seconds.
Move from down to up and laterally.

*When Galvanic pin working, every 5 seconds with "BI", mainly remind nursing to move point by point. Every 10 minutes the alarm, a beeping sound will get off and the appliance will be turned off automatically. Then you can begin with another parts of your body by pressing the ON/OFF button on. If you would like to stop the machine before that, please press ON/OFF button any time to turn off machine.

The gel-like essence is recommended for this is helpful for ultrasound wave transmission.

6

Daily caring

Do not use the machine more than 3 times in one week, and do not work more than 30 minutes each time.

Step 1:

Deep cleaning skin procedure

- Thoroughly remove make-up and impurities first.
- Apply appropriate portion cleaning cream on your face and the probe surface.
- Switch ON, with ionic- function, and select the intensity you want.
- According to the diagram at right for Massage, Slowly moving the probe that along indicated route from jaw to the root of your ears with 5 seconds, other parts of are the same way.
- Clean face with water and clean machine.



Step 2:

Induct Nutrition and shape your face


- Apply appropriate portion essence on your face and the probe surface.
- Switch ON, with ionic- and ultrasonic function, select a intensity make you comfortable.
- According to the in page 8 for Massage, Slowly moving the probe along indicated route from jaw to the root of your ears with 5 seconds, other parts of are the same way.

Step 3: Photon caring


- Switch ON, Select "Photon" Function, and choose one color of photon you want:
Red: Promote collagen protein growth. Anti-wrinkle
Blue: Sterilization and acne removing
Green: comfortable, make skin relaxed
- After photon selection, Using photon irradiated massage in need of treatment sites with point-by-point moving.

Every treatment place, need about 10-30 minutes.


7




1. Move from down to up and geometrically.




2. Slowly moving the probe from jaw to the root of your ears.




3. From mouth corner to front of ear.



4. From nose position to brow corner.



5. Nose position.



6. From nose to temple.

8

Photon ultrasonic anti-wrinkle machine

Remove Wrinkle Nursing

(For wrinkles at canthus, corner of mouth and forehead.)

Step 1: Deep cleaning skin

(The same as Daily caring Step 1)

Step 2: Anti-wrinkle

- Apply wrinkle care essence on the wrinkle skin and the probe surface, make them moist enough.
- Setting Ultrasonic function, and switch on the red light.
- Choose a comfortable intensity from "LEVEL" button, and Refer to the in page 8 Anti-wrinkle massage;

(10 minutes power-off automatically, you can open it again, if you want to continue doing treatment)

Dilute the spot, even skin color, whitening

Step 1: Deep cleaning skin

(The same as Daily caring Step 1)

Step 2: Select Ion- Function

- Induct-in Concentrate on spot area, Apply appropriate portion essence, such as VC, on your face and The probe surface. During the massage, the probe should staying 20-30 seconds (5-6 Tones) in each position;
- Switch ON, with ionic- and ultrasonic function, then select the intensity you want.
- At the same time can choose the blue light massage according to page 8, with photon functional effect is more obvious.

9

