

3-in-1 Multifunctional Eye-care Reviving

User Manual

- Thank you for purchasing this 3-in-1 Multifunctional Eye-care Reviving
- Please use under guidelines in this book and be well aware of the safety information listed. Keep this book for future reference
- This device is designed to provide routine eyes and skin cares. Any other ways of use and sale in connection with medical purpose are not supported

3-in-1 Multifunctional Eye-care Reviving

It is a small pen-shape beauty device that provides a whole set of eye refining treatments. It has a rolling ball on one end and ultrasonic steel head on the other. Through only two steps, you can use it any time in any places to easily:

1. Replenish nutrition into skin for deep moisture and hydrating
2. Firm and enhance skins
3. Repair crows feet and fine lines
4. Lighten dark circles, eye pouches and puffiness
5. Cleanse and detoxify the lymph system

CONTENTS

1.Product Details	01
2.Safety Instructions	02
3.Parts Introduction.....	03
4.Functions & Operation Instructions	04
5.Before Use	05
6.User Guide	06
7.After Use	08

Product Details

Name: 3-in-1 Multifunctional Eye-care Reviving
Input power: AC 100V – 240V
Output Power: DC 5V 1,000 mA
Dimensions (cm): 16 x 4.8 x 3.5
Weight (g): 118
Diameter of Ultrasound Head (mm): 25
Timer setting: Automatic shut off after 5 minutes
Ultrasonic Frequency: 3M Hz ± 5%


Safety Instructions

- This device is designed to provide routine eyes and skin cares. Use for any other purpose is not supported.
- Please read this chapter carefully before use. Follow the guidelines while use and control the frequency of use under instruction of this manual. Please shorten the operating time for the first time of your use if your skin is very sensitive.
- Please stop using if you feel any uncomfortable. Consult your doctor right away.

Do not use this device

- With any medical machines such as Cardiac Pacemaker, Oxygenator, electrocardiograph, etc.
- In any of circumstances below:
 - During pregnancy or menses
 - Have allergic dermatitis or hemorrhage
 - On Irritating skin or redness caused by sunburn
 - On skin ulcers or inflammations
 - On cosmetic operated body parts
 - Have telangiectasia or disordered liver function caused by long term uptake of steroid and hormone


⚠ Operation & Maintenance Precautions

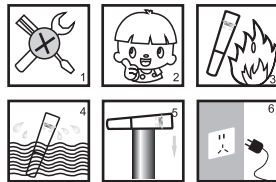

NOTE

- Do not disassemble it by yourself
- Keep it away from children
- Keep it away from fire
- Keep it away from humidity and water
- Do not drop or collide it
- Unplug the electricity cord after use

Parts Introduction

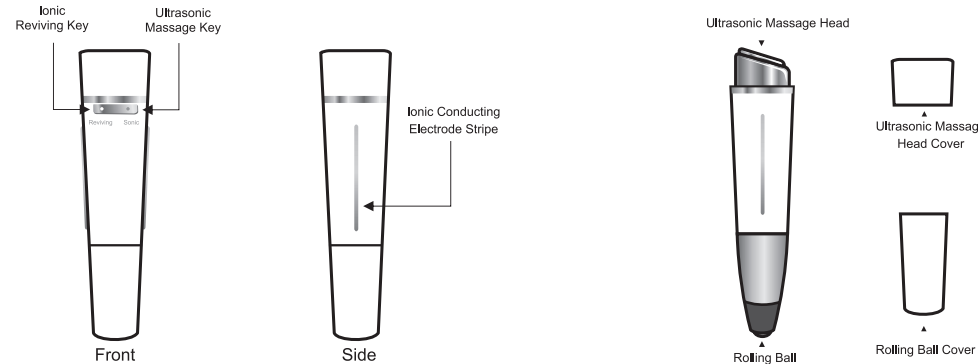

Functions:

The technologies we choose include 3M Hz ultrasound, ion conducting & heat massage, which are clinically proven to be safe and effective, and fit all skins. Just several minutes, awake youth and beauty of your eyes from inside.

- 3M Hz ultrasound: cleanse and detoxify the lymph system; take away surplus water and fat around eyes. Firm and lift your skin with visible results after use.

- Ionic conducting treatment: replenish your eye care product deep into skin to enhance its effects; moisture and hydrate your skin; lighten dark circles and eye pouches.

- Heating massage: the built-in $40^{\circ}\text{C} \pm 2^{\circ}\text{C}$ warming system can relieve eye strain, speed up blood circulation of capillaries around eyes, thus reduce fine lines, crows feet and make your eyes shine.

- Acupoints stimulating: while moving around your eyes, it can stimulate the acupoints around, refresh your skin and strengthen diminishing of puffiness and pouches.

Operation Instructions:

- **ON/OFF:** long press on 'Reviving' key or 'Sonic' key for power on/off.

- **Select Function:**

Ionic Conducting with Heating Massage: 'Reviving'
Ultrasonic Treatment: 'Sonic'



Confirmation of Ultrasound:

Lay the ultrasonic head flat with head upwards, drop a water onto the head, and press 'Sonic'. The water drop will vibrate when the ultrasound is on.



Confirmation of Ionic Conducting with Heating Massage:

Heating Massage:

Press 'Reviving', and the rolling ball will slowly raise its temperature to around 40°C in 3 minutes (28°C indoor temperature for reference). Touch the rolling ball onto skin and the heat can be felt.

Ionic Conducting:

You may feel tiny tingle or numbness when ionic conducting is on. Different skincare products and different skin types can lead to different sensations.

- Shut Off: both 'Reviving' and 'Sonic' are set 5 minutes for automatic shut-off. The device will shut off after a beep.
- Charging: please charge before your first use and every time the battery runs out of electricity.

Charge at least 3 hours before first use to activate lithium battery inside.

In charge: sync flash of two signal lights

Charge complete: continuous lighting of two signal lights

Out of electricity: circulation flash of two signal lights, automatic shut off after 10 times of flash circles.

Charging Time	Overall Operation Time
About 3 hours	About 2 hours

(The above operation time is based on immediate use after full charge)

Before Use

Frequency of use:

Mode	Frequency of use
Reviving	2-3 times/week, 3-5 minutes/time
Sonic	2-3 times/week, 5-10 minutes/time



Quick User Guide:

Visible results soon after use to take care of all skin concerns around eyes (dark circles, pouches, puffiness, crows feet, etc.), and can also be used on face, helping to regain elasticity of your skin.

Step 1:

Remove makeup and clean the skin

Step 2:

Apply your eye care product evenly around eyes, and conduct the nutrition deep into skin through REVIVING function (see next pages for more details) to diminish black circles and pouches.

Step 3:

Apply gel evenly onto face and do massages with SONIC function (see next pages for more details) to take away extra water and fat, reduce fine lines and crows feet, as well as clean and detoxify the lymph system.

Tips:

1. Place the ultrasonic head softly onto your skin, and ensure it contacts your skin tight during use.
2. If the rolling ball or ultrasonic head is not lubricated enough to move smoothly on your skin, apply more essence or gel onto it.
3. Keep it away from eyeballs, and avoid long time stay on temples

User Guide for REVIVING function:

1. Cleanse the skin.



2. Apply your eye care product around eyes, following the direction from the inner eyes corner to the temple.



3. Unplug the cover on the end with rolling ball, and press 'Reviving' button.



4. Hold the device with your hand touching the electrode stripes on both sides.
(Moist your hand with wet tissues first if it is dry.)




5. Slowly move the device back and forth from the middle of your lower eyelids to the inner and outer corners of your eyes, replenishing nutrition deep into skin.



User Guide for SONIC function:

1. Unplug the cover on the other end with ultrasonic head.

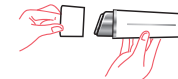

2. Apply your skin care product / gel evenly onto skin as well as the ultrasonic head.



3. Start 'Sonic' function, press the ultrasonic head softly on skin and massage following paces below:



Push the device slowly upwards along the curve of chin to behind the ears.



Push the device slowly downwards along the neck to the shoulder to cleanse and detoxify the lymph system and promote the body circulation.




Push the device slowly upwards following arrow directions in picture.
(Repeat each step for 5-10 times.)



After use:

- Take off the protection cover around the rolling ball by rotate the cover counterclockwise, clean the cover and the rolling ball with water flow, drain them and rotate the cover back.
(Keep the other parts of device away from water while cleaning)

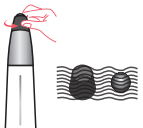

- Swipe off remained skincare product / gel on the ultrasonic head with wet tissue or rags.



- Keep the device clean and dry after use.