



Dedicate to easying quick effective beautify skin in your home



SKIN CARE MACHINE

Operating Instructions

Ionic-Ultrasonic Skin Care Machine

- ◆ Thank you for purchasing this Ionic-Ultrasonic Skin Care Machine
- ◆ Please read this Operating Instructions manual carefully before use, and keep this book in a safe place for future reference
- ◆ This product is designed to take care of normal skin and to beautify it. Any other usage and sales in connection with medical purpose are not supported

Contents

- Specification
- Safety precautions
- Important notice
- Skin sensitive test
- Parts Introduction
- Operating notice
- Preparation before using for the first time
- Operating Methods
 - Ionic cleansing caring
 - Ultrasound Facial Massage
 - Suggested operation process
- Emphasized treatment procedure around eye

Specification

Product type: ZHF-CM-3
Battery: DC 6V/12V
Current: Max. 20 mA
Adaptor input: AC 100-240V
Adaptor output: DC 15V/ 0.4A
Power consumption: Max. 300 mA
Dimension: 95mm×55mm×45mm
Ultrasonic frequency: 3MHz ±5%
Weight: 49g

Safety precautions



Note: Please read following information carefully before using your machine. This product is well-designed precision instrument. In order to have your device functioning perfectly, please read the below information thoroughly and follow all the instructions. This will help you to ensure the warranty validity and to extend the life cycle of your product.

- ◆ 1. Store your product always from the humidity and do not use it under water. Otherwise the electronic components of your product be corroded or have short circuit.
- ◆ 2. Never attempt to disassemble the product yourself. All the product components can only be maintained and repaired by professional technicians.
- ◆ 3. Don't use any adaptor other than the one included in the package. Using incompatible adaptor may result in malfunction and will void the warranty validity.
- ◆ 4. Any people with the following conditions must avoid using this product: patients with malignant, women during her pregnancy or period, patients with dermatosis of face, on thrombosis, urgent patient with unidentified diseases, and children and elders. Forbid to use at special position, such as ear, throat, eyeball.
- ◆ 5. Keep this product away from children. Children shall not play with this product as the disassembled components may be swallowed by the children accidentally.




Important notice


- ◆ To ensure that this machine functions well, you can drop a tiny amount of water on the ultrasonic wave generator surface to see the water bobs. If the water bobs, it means that this machine is under normal function.
- ◆ This machine should be used under normal conditions. Remember to turn off the power and unplug the cable when you stop using skin. This will avoid over-heating problem caused by the air vibration and to prevent any possible burning of shortening product life cycle.
- ◆ Don't use the machine without putting any skin care product on your skin-this will cause wave seriously reflection and result in discomfort or burning of your skin. Gelatinous skin care product is recommended.
- ◆ Avoid using this machine in one area for too long. Move it around on your skin to avoid focusing effect caused by the wave energy. This effect may result in discomfort or burning of your skin.
- ◆ When you are using this machine, always adhering the ultrasonic wave generator completely to your skin. An incomplete contact will cause a big difference of impedance and disperse energy on the interface. This will also cause an over-heating problem and result in a discomfort or burning of your skin.
- ◆ When you finish using the machine, wipe off the medium left on the ultrasonic wave generator. A soft cloth is recommended to prevent any possible scratch.
- ◆ Do not use water and any other chemical, cleanser or washing-up liquid to clean this product. This will cause a malfunction or shortened product life cycle.

Please do allergy test before using for the first time.(as follows)


- ◆ Cleanse your wrist.
- ◆ Apply a portion of the special essence on the cleansed area, let it dry naturally (During the test, please don't touch water nor wash the essence away). If you experience allergy or any skin abnormalities, stop the test immediately and contact your doctor or nearest hospital straight away.
- ◆ If after 48 hours, no allergy or any other forms of skin abnormalities appear, you may enjoy using this special unique essence and equipment regularly.

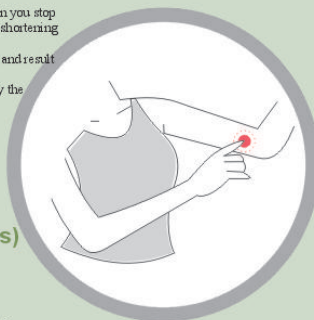

Parts introduction

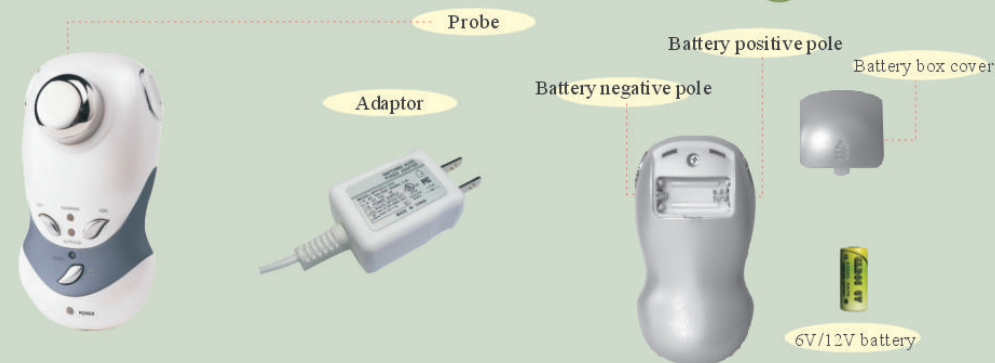

Operating Panel Instruction

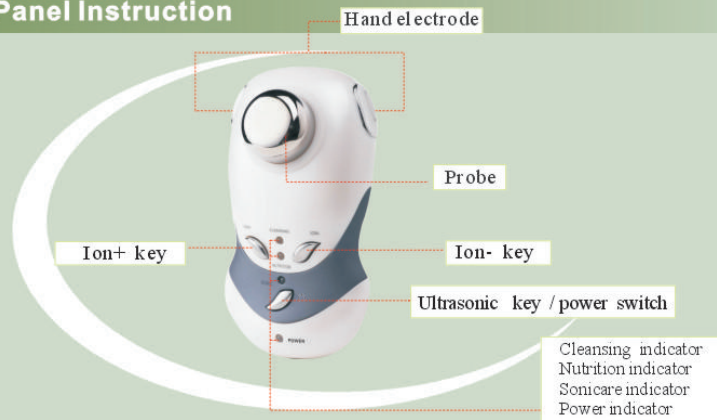

Operating instruction:



1. For your convenience, this machine provides two kind of power sources : battery and adaptor.
 - (1). For portable, this machine could be powered by battery which could only act at "Ion-" and "Ion+" .
* When the battery voltage is low, the four indicators will flash orderly (after 10 second the power will turn off automatically). Please replace all the batteries with new ones.
 - (2). When the machine using adaptor ,the machine will automatically use adaptor power no matter with or without battery and all functions, such as Ion-, Ion+, sonic, sonic++ionic, sonic, sonic-tionic can working.
2. All function have two level low and high There is one "Bi" sound for low level and one "BiBi" sound for high level.
3. There is one sound for every five seconds, "Bi" for low level and "BiBi" for high level.
4. After 10 minutes ,power will be off automatically.
5. Press power button two seconds can switch off machine at any time and any state.

Operating method

Before operating

1. power confirmation

(1) Remove the battery box cover, install 1pc 6v or 1pc 12v battery. Align the poles (+ and -) properly when inserting the batteries.

(2) Operation Procedure for testing ultrasonic working

- ◆ Connecting adaptor, place the appliance on a flat surface and put several drops of water on the ultrasonic waves generator.
- ◆ Press power switch on. Press ultrasonic key one time (ultrasonic low level) or two times (ultrasonic high level).
- ◆ If the water on the ultrasonic waves generator vibrates, the ultrasonic waves are produced.



▲ Connecting adaptor, place the appliance on a flat surface and put several drops of water on the ultrasonic waves generator.



▲ If the water on the ultrasonic waves generator vibrates, the ultrasonic waves are produced.

2. Remove any jewelry i.e. earrings, rings, necklace, watch and so on.
3. Cleanse your face and remove any make-up gently and thoroughly with cleansing lotion or other facial cleansing liquid.



Operating method



I. Ionic cleansing (Ion +) caring

1. Apply appropriate portion essence on your face and avoid eye area. Please ensure the treatment area is well covered with the essence.

(ensure your face and the cosmetic cotton tissue are thoroughly moisturized) :

2. Turn on power (press power switch one time) :

3. Operation order:

- (1) Push Ion+ button once --the machine works at Ion+ of low level meanwhile the cleansing indicator on.
 - (2) Push Ion+ button once more and the machine works at Ion+ of high level.
 - (3) Push Ion+ button third times and Ion+ turn off.
 - (4) In any level of Ionic out, you can turn to Ion- function by pushing the Ion- button.
 - (5) In any function (Ion+, or Ion-) and in any level, you can turn on sonic function too, then sonic care indicator on and ionic function change to low level automatically.
4. When you use Ion+ or Ion- function, you must hold the Hand Electrode during operating. Glide it slowly in the order indicated in the diagram.
5. After 5 minutes caring time, power will be off automatically.



When you use Ion+ or Ion- function you should hold the Hand Electrode with thumb and second finger, otherwise the Ion- and Ion+ will not be acted.

Operating method

II、Nutritious caring

It should be done after Ion+ caring.

1、Apply appropriate portion essence on your face and avoid eye area. Please ensure the treatment area is well covered with the essence. (ensure your face and the cosmetic cotton tissue are thoroughly moisturized) ;

2、Press power switch one time ;

3、Operation order

(1) Push Ion- button once --the machine works at Ion- of low level meanwhile the nutrition indicator on.

(2) Push Ion- button once more and the machine works at Ion- of high level.

(3) Push Ion- button third times and the Ion- turn off

(4) In any level of Ion- You can turn to Ion+ function by pushing the Ion+ button.

(5) In any function (Ion+, or Ion-) and in any level ,you can turn on sonic function too and Sonicare indicator on and ionic function change to low level automatically.

4、When you use Ion+, or Ion- function, you must hold the Hand Electrode during operating. Glide it slowly in the order indicated in the diagram.



working upwards motion.
Caring both right and left parts symmetrically.
Glide slowly (It should take 5 seconds from chin to the bottom of ear) Diagram 1

Operating method



III、Ultrasonic caring operating

Basic method

1、Apply appropriate portion essence on your face and avoid eye area. Please ensure the treatment area is well covered with the essence, and apply appropriate portion essence on the probe.

2、Connect the machine with adaptor. Press power button one time ,power indicator on (strongly emphasize, only adaptor can drive ultrasonic function)

3、Operation order

(1) Push Ultrasonic key once --the machine works at sonic of low level meanwhile the sonic indicator on .

(2) Push Ultrasonic key once more and the machine works at sonic of high level.

(3) Push Ultrasonic key third times and sonic function turn off.

(4) In any level of sonic ,you can turn on Ion- or Ion+ function by pushing the related button, and related indicator on. Glide it slowly in the order indicated in the diagram. Hold the Hand Electrode during ionic care.

4、Wash off remaining beauty gel on face after usage.

Operating suggests

I、Basic operating: twice a week in the beginning and once a week after two months.

- 1、 Ion+ and ultrasonic function for deep cleaning (3-5min)
- 2、 Then Ion- and ultrasonic function for nutrition (3-5min)
- 3、 At last Ultrasonic operating Enhance absorption of nutrition (3-5min)

II、Treat macula by the help of whitening, moisture keeping, wrinkles preventing essence.

- 1、 Ion+ function for deep cleaning(3-5min)
- 2、 Then Ion- and ultrasonic function for nutrition induct-in (3-5min)

Twice/ third times a week in the beginning , then decrease care time gradually and distinct effect appears.

Strongly emphasize: (do not use this appliance directly on the eyebal)




Every woman has had one of those mornings when the mirror simply doesn't hide the truth. Lack of sleep, illness, irritation , stress or too much Caffeine all take their toll on your skin. For many, this translates into puffy eyes or black or crow's feet.

The skin around eye is much more thin than other skin , it need more superficial and careful care. There is 3 million vibrations per second on the probe area and 3 million vibrations per second can not reach deep skin.

There are mainly three type of problems for eye surrounding:

black eye, puffy eye , and crow's feet .

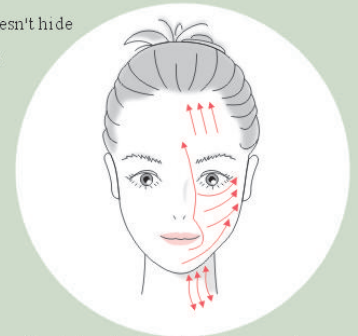

Should not stay too long at forehead and ear area

Emphasized treatment procedure around eye.

I、Operation of crow's-feet :

The temple can not be operating for a long time .

Saggy skin caused by lessen collagen for the normal aging process.

So you should care it carefully.

Caring procedure:

- 1、 Ion+ function for deep cleaning(1-3min)
- 2、 Ion- and ultrasonic function for nutrition (2-4min) (the beauty gel should contain Hyasol or phytokine to moisturizing skin or anti-age)
- 3、 At last gently apply a eye mask to contain the effect .



II、 Puffy eye---- forbid massage the skin upon eyeball.

Puffy eye ,also known as “bags,” can be caused by poor circulation. To speed the blood circulation and balance metabolism, so the beauty gel should have Dipeptide-2, Hyasol and so on.

Caring procedure:

1、 Ion+ function for deep cleaning(1-3min)

2、 Ultrasonic operating:

- Along indicated route from jaw to the root of your ears.
- Along indicated route from root of your ears to supraclavicular fossa.
- Along indicated route from supraclavicular fossa to armpit
- Along indicated route from shadow under eye to the temple (should not stay too long on temple area).
- Along indicated route from the temple to opisthotic .

