




SR-5

Salon Quality Skin Technology



BIO-SCIENCE™
RESEARCH

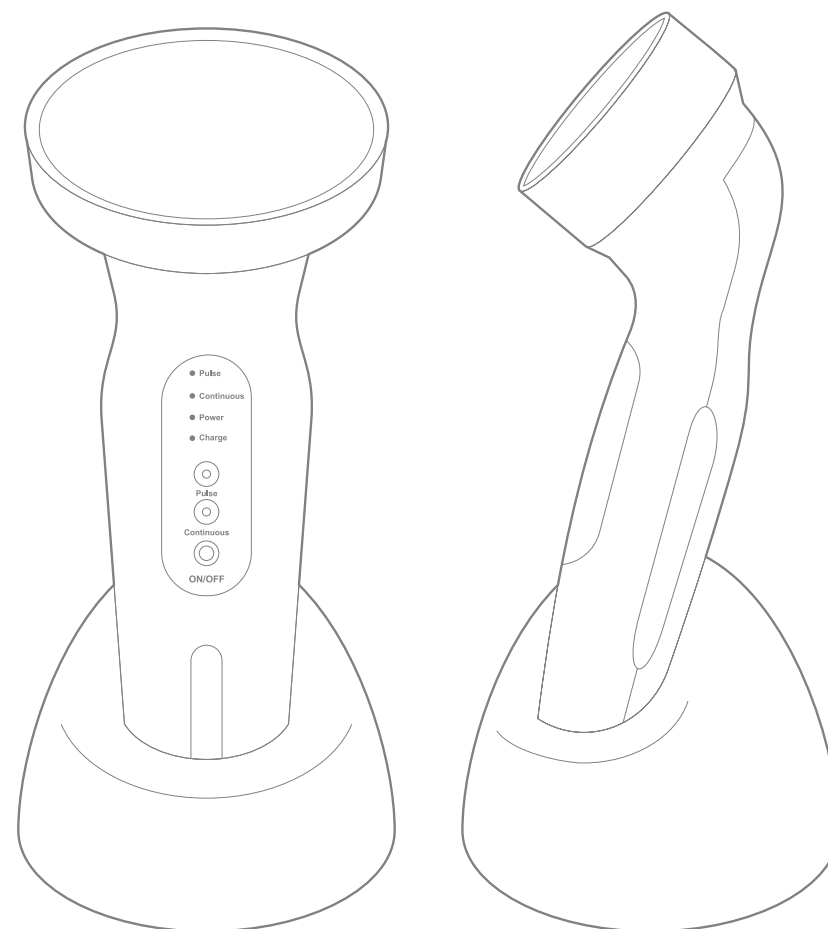



Congratulations!

Specification


Product name: SR-5 Skin Technology

Adaptor input power: AC 100-240V

Adaptor output power: 5V \pm 5%, 1A

Adaptor frequency: 50/60 Hz

Photon wavelength: 620-625 nm

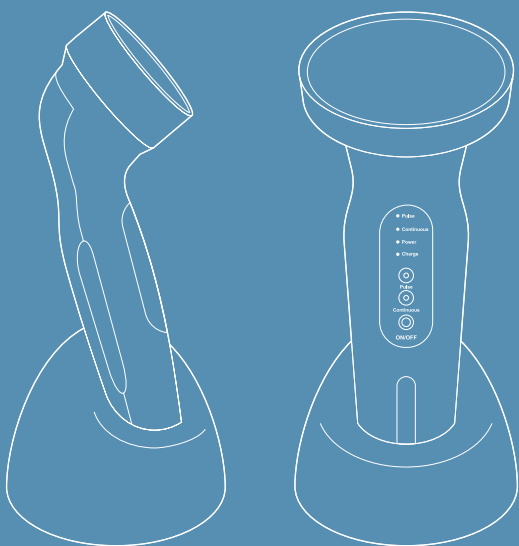

Yes, now you can get rid of wrinkled, sagging skin and look many years younger! Visibly younger!

Without painful injections, damaging laser therapy or expensive trips to the doctor. At last, a clinically proven, safe and painless solution gives you salon professional results... in the comfort of your own home.

Plastic surgeons around the world are now using light therapy - called photo-rejuvenation - to successfully reduce wrinkles and fine lines. Painful treatments like Botox® paralyse the facial muscles, creating only the "effect" of softer skin. And in many cases giving your face that bloated, plastic doll-like appearance. Yuck!

But clinical studies now show that photo-rejuvenation treatments are proven to increase your skins collagen and elastin levels - naturally removing wrinkles from the inside out.

Clinically proven to be pain free and safe for all skin types, these specially calibrated light waves (620nm) penetrate the top layer of the skin causing the cells to produce collagen and elastin. And more collagen and elastin means your skin naturally firms and tightens... visibly reducing the appearance of fine lines and wrinkles. Leaving your face feeling smoother and more vibrant.

Photo-rejuvenation treatment has also been shown to reduce age spots and freckles, visibly lower skin redness and irritations and to accelerate skin repair.

Safety Precautions

Do not use if you:

- * Have had plastic surgery, fracture repair surgery, metal, silica-gel or any kind of plastic implanted in your face
- * Have inflamed skin due to over-exposure to sunlight
- * Have dermatitis caused by cosmetics
- * Have itching or hot spots
- * Have spots caused by bites

Operation and Maintenance Precautions

1. Store product away from humidity and do not use under water, otherwise the electronic components of your product may be corroded or short circuit.
2. Never attempt to disassemble the product yourself; product components can only be maintained and repaired by professional technicians.
3. Don't use any adaptor other than the one included in the package. Using incompatible adaptors may result in malfunction and will void the warranty validity.
4. Keep product away from children. Children should not play with this product as the disassembled components may be swallowed accidentally.

Note: Please read this information carefully before using your machine. The product is a well designed precision instrument. In order to have the device functioning perfectly, please read the information thoroughly and follow all the instructions. This will help to ensure the warranty validity and to extend the life of the product.

2



Important safety notice



Do not open or disassemble the device yourself.



Keep product away from children




Keep away from open flames.



Do not use in or near water. Please store in a dry place.



Dropping the device may cause it to stop working.

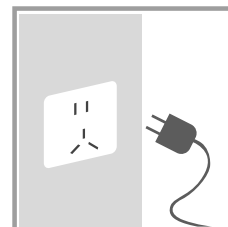


Please un-plug docking base after use.

3



Device

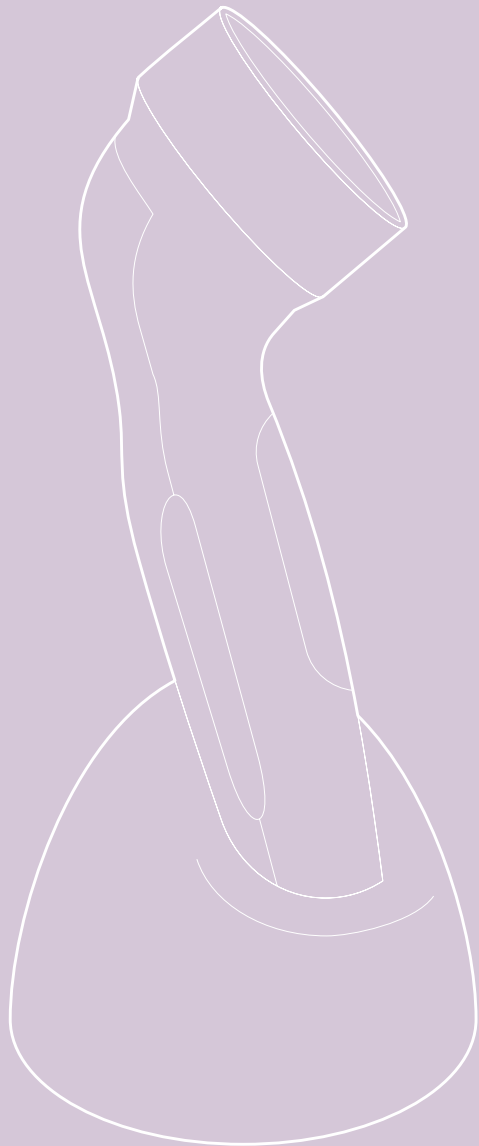

4

Overview

Salon Quality
Skin Technology




LED Array Contact Surface


Pulse Mode Indicator
Continuous Mode Indicator
Power Indicator
Charge Indicator

Pulse Mode Selector
Continuous Mode Selector
Power ON / OFF Selector

Recharging Docking Station

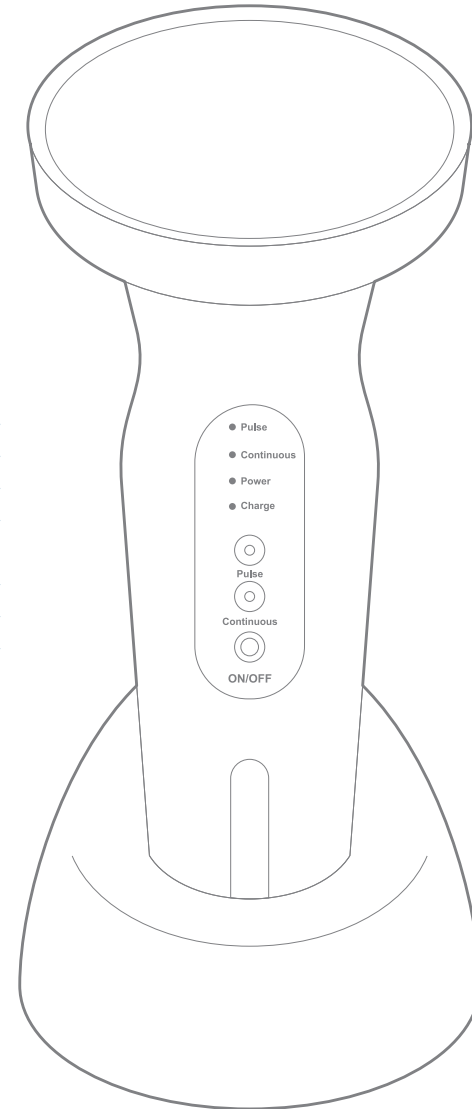

5

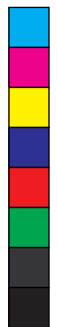



Getting ready for use:

Using around your eyes


While the LED lights are safe to use on your face, it's recommended that you don't point the light directly into your eyes. When using the device, please use the goggles provided, or close the eye closest to the device when using.

Device Care

To clean the device, simply use a small amount of rubbing alcohol and a cloth.

Charging


To charge the device prior to use, simply plug the docking station base into the power supply and place the device into the docking station base.

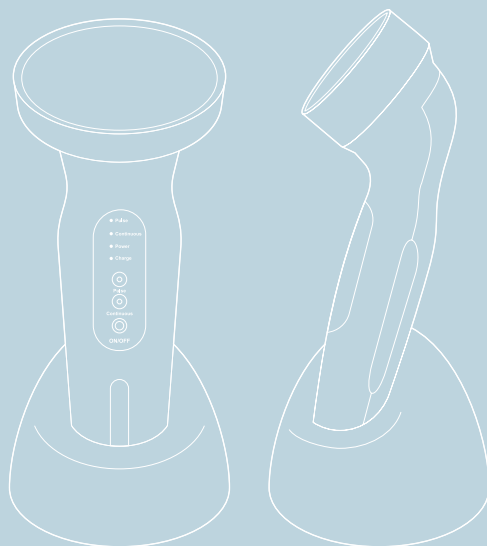

6


TIPS FOR USING THE DEVICE

When using the device, press the round Contact Surface gently onto the skin. The Contact Surface should be flat on the skin, not on an angle.

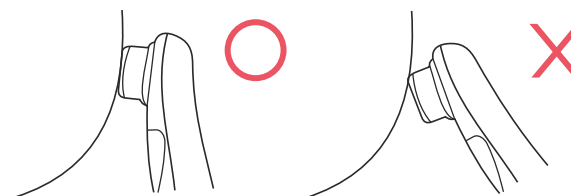

Before You Start... It's Photo Time!

Before you use SR-5 Skin Technology for the first time, take a photo of your face up close so you can see your progress. In order to get an accurate 'before' and 'after' photo, make sure you take the photos in a well lit room each time. Every 8 weeks, take another photo in the same room and lighting to get an accurate record of your progress.

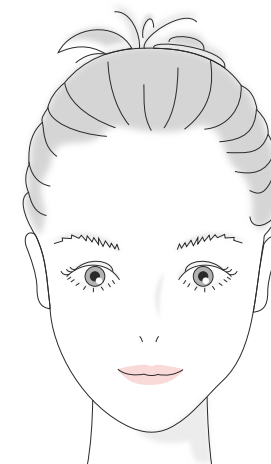

7

Treatments

QUICK START GUIDE

Photo-rejuvenation treatments are so gentle and relaxing that many people enjoy treating their wrinkles daily while watching TV, chatting on the phone or even reading. It's really that simple and easy!

While many people report visible improvements in just a few weeks, it is recommended to continue use for at least 8 weeks for maximum effect.

Step 1: Cleanse Your Skin

Cleanse your face thoroughly using your usual skin care products.

Step 2: Skin Activation (Continuous Mode)

Turn on the device and select "Constant" Mode. This provides a constant output of light energy that lasts for 3 minutes.

Cover your eyes with the goggles provided.

Starting with your forehead, touch the device gently to the skin. The light will automatically activate as it touches the skin.

Following the treatment diagrams, slowly move the device around the entire face.

Spend 30 seconds on each of the 6 treatment areas. Use small circular motions over each treatment area to ensure enough light penetrates the skin surface and activates the cells.

The entire Skin Activation treatment will last 3 minutes before the device automatically turns off.

Step 3: Skin Treatment (Pulse Mode)


Turn the device on and select "Pulse" Mode. This provides short specific pulses of light energy that are shown to enhance treatments.

Cover your eyes with the goggles provided.

Starting with your forehead again, touch the device gently to the skin.

Following the treatment diagrams, hold the device in place over the desired treatment area for 4 minutes. The device will automatically switch off after 4 minutes.

Repeat this process until you have treated all the desired areas.

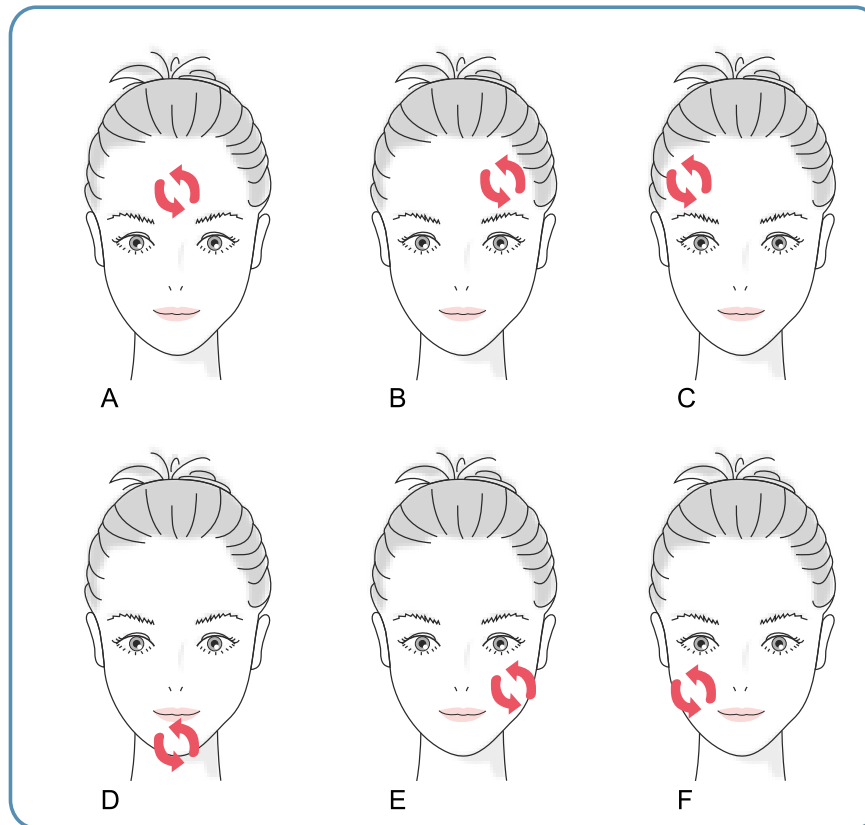

Step 4: Other body parts

You can also use the device on sagging skin and fine lines around your neck, chest and other body parts.

