

[Guarantee Certificate]

Product model:_____	Purchase date:_____
Purchase place:_____	Sales Dept. _____

One year Repair Guarantee

In the unlikely event that your product develops a fault within the first one years of purchase:
(Please return your product along with your receipt and the original guarantee certificate and we will repair it free of charge provided you have operated it according to the Operating Instructions book)

What is not covered?

1. Breakdown or damages due to
- (1)Not following the operating instructions.
 - (2)Incorrect or unauthorized repair
 - (3)Accidental damages during your own transportation, or you drop it on the floor or you deliberately damage or neglect the product after purchase.
 - (4)Fire, earthquake, weather such as flood, lightning, high wind.
2. Please keep your guarantee certificate safe, as no duplicated certificate will be issued.

PHOTON- ULTRASONIC WAVES SKIN CARE MACHINE

PHOTON- ULTRASONIC WAVES SKIN CARE MACHINE
Operating Instructions

- Thank you for purchasing this Photon- Ultrasonic Waves Skin Care Machine.
- Please read this Operating Instructions manual carefully before use, and keep this book in a safe place for future reference.
- This product is designed to take care of normal skin and to beautify it. Any other usage and sales in connection with medical purpose are not supported.



Table of contents

- Specification
- About charging
- Safety Precautions
- Important notice
- Parts introduction
- Quick user guide
- Photon treatment
- Warranty certificate

Specification

- Photon-Ultrasonic Waves Skin Care Machine
- Adaptor input: AC 100-240V
- Adaptor output: DC 5V 1000mA
- Dimension : 17.5*4.7*4.4cm
- Weight : 150g
- Ultrasonic frequency: 3 MHz \pm 5%

About charging

The capacity of battery in storage and transportation is only 30%, so please you charge it before usage. The normal charge time is 3 hours.



- Recharge with your computer USB source by connecting machine to computer with USB line.

- Recharge with adaptor by inserting adaptor plug into wall socket.

- Recharge with charge socket with computer or adaptor.

Operation Cautions



- The ultrasonic waves generator must be contacted with the skin.
(You MUST press the ultrasonic waves generator well against the skin and ensure there is no gap between them.)



- How to glide the ultrasonic waves generator:

1. Glide slowly it takes 5 seconds from jaw to the bottom of your ear.
2. Glide upward motion, on both right and left areas symmetrically.



It normally requires 4-6 seconds each time.

- Please apply enough essence on the surface of ultrasonic waves generator to ensure its smooth operation. If requires please add more essence on.

photon treatment

photon caring: Daily caring .

Select red light, nursing face , especially wrinkle area.
move detector every 5 seconds.

Recommended photon caring time is 10-20 minutes

photon caring: Dilute the spot , even skin color, whitening skin

Select blue light, nursing face , especially spot area.
move detector every 5 seconds.


Recommended photon caring time is 10-20 minutes

photon caring: Acne (one time every night)

Lymphatic detoxing

1. Select ultrasonic function and intensity, perform Lymphatic detoxing according to next route.

Use blue light to care acne area 15-20 minutes.
Use green light to care 10 minutes.

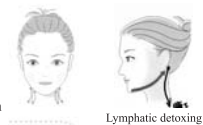

Lymphatic detoxing



Massage around acne.

Safety precaution

Note: Please read following information carefully before using your machine

This product is well-designed precision instrument. In order to have you device functioning perfectly, please read the below information thoroughly and follow all the instructions. This will help you to ensure the warranty validity and to extend the life cycle of your product.

- 1.Store your product always from the humidity and do not use it under water. Otherwise the electronic components of your product be corroded or have short circuit.
- 2.Never attempt to disassemble the product yourself. All the product components can only be maintained and repaired by professional technicians.
- 3.Don't use any adaptor other than the one included in the package. Using incompatible adaptor may result in malfunction and will avoid the warranty validity.
- 4.Any people with the following conditions must avoid using this product: patients with malignant, women during her pregnancy or period , patients with dermatosis of foci on themorrhage, urgent patient with unidentified diseases, and children and elders.
- 5.Keep this product away from children . Children shall not play with this product as the disassembled components may be swallowed by the children accidentally.

Important notice

- 1.To ensure that this machine functions well, you can drop a tiny amount of water on the ultrasonic wave generator surface to see the water bobs. If the water bobs, it means that this machine is under normal function.
- 2.This machine should be used under normal conditions. Remember to turn off the power and unplug the cable when you stop caring skin. This will avoid over-heating problem caused by the air vibration and to prevent any possible burning of shortening product life cycle.
- 3.Don't use the machine without putting any skin care product on you skin-this will cause wave seriously reflection and result in discomfort or burning of your skin. Gelatinous skin care product is recommended.
- 4.Avoid using this machine in one area for too long. Move it around on your skin to avoid focusing effect caused by the wave energy. This effect may result in discomfort or burning of your skin.
- 5.When you are using this machine, always adhering the ultrasonic wave generator completely to your skin. An incomplete contact will cause a big difference of impedance and disperse energy on the interface. This will also cause an over-heating problem and result in a discomfort or burning of your skin.
- 6.When you finish using the machine, wipe off the medium left on the ultrasonic wave generator. A soft cloth is recommended to prevent any possible scratch.
- 7.Do not use water and any other chemical, cleanser or washing-up liquid to clean this product. This will cause a malfunction or shortened product life cycle.

Please do allergy test before using for the first time.

The Simple Allergy Test:

- Cleanse your wrist.
- Apply a portion of the special essence on the cleansed area, let it dry naturally (During the test, please don't touch water nor wash the essence away). If you experience allergy or any skin abnormalities, stop the test immediately and contact your doctor or nearest hospital straight away.
- If after 48 hours, no allergy or any other forms of skin abnormalities appear, you may enjoy using this special unique essence and equipment regularly.



-2-

Operation to tighten face

- Operate from the corner of mouth towards the area near eyes.




Operation to pull up the end of mouth

- Operate from the corner of mouth in upward motion



Eczema and acne care

- Do not operate on the area where there are eczema, macula, acne, sunburn etc directly. However, you may care the surrounding skin area to aid in the natural healing of skin disorders.

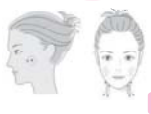

Operation to tighten double chin

- Glide it from chin to bottom of ear, from Chin to the corner of your mouth. And always
- Glide it from lower parts to the upper parts.



About T position caring

- As indicated in the diagram, Glide your machine in nosewing and on the bridge of your nose for 1-2 minutes.



Dispels the black eye socket, eye bag

- By in corner of the eye underneath to temples
- Along face edge to root of the ear behind
- From root of the ear behind to collar bone spot
- From collar bone spot to armpit
- It's better for taking 20 times



-7-


Special method to care your face

(Operation time: 2 minutes)

The gravity of our earth, increasing age, getting fat, all these factors are causing your skin loosening, the prolapsed skin and having a double chin. For the following operations, it takes 10-20 for each area.


Please apply slight strength on the ultrasonic waves generator to lift up the skin and muscle following the arrow direction indicated in the diagram. After this basic skin-care method, now you may use special skin care methods to operate.

Red light: excitement Green light: comfortable Blue light: composure

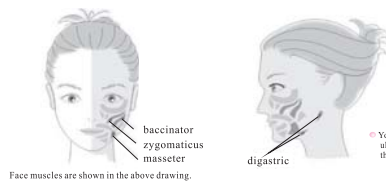

- 6 -

Parts introduction

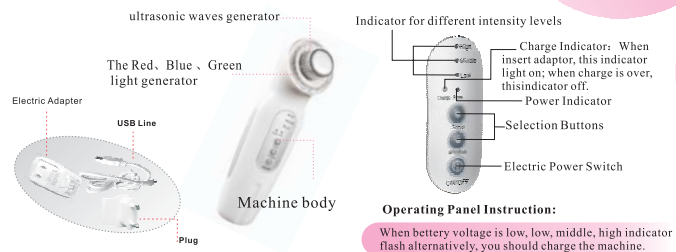

Operating Panel Instruction:

When battery voltage is low, low, middle, high indicator flash alternatively, you should charge the machine.

- 3 -

QUICK USER GUIDE

Method to care your face

- 1. Thoroughly remove make-up and impurities first. The ultrasonic waves skin-care treatment must be performed on the cleansed skin.
- 2. Turn power on and then the indicator light is on.

- 3. Apply appropriate portion essence on your face and avoid eye area. Please ensure the treatment area is well covered with the essence. (Note: We highly recommend you to use our designated essence for its best performance on conducting ultrasonic waves)
- 4. Apply appropriate portion essence on the ultrasonic waves generator.



Suggestion: You'd better use the designated essence to transfer ultrasonic waves.

-4-

Basic method to care your face



Step 1 it requires 8 minutes.

Using ultrasonic waves generator to massage skin equally.

It requires 4 minutes. First, Using ultrasonic waves generator to care the right side of your face for 4 minutes (Please use a timer or your watch as a timer). Then working on the left side of your face for 4 minutes.

Please operate the appliance in the order indicated in the diagram, working upward motion (about three times).



1. Along indicated route from jaw to the root of your ears.



2. Along indicated route from jaw to the front of your ears.



3. Along indicated route from the corner of mouth to the top of your ear.



4. Along indicated route from the bottom of nose towards the area near your eyebrow.



6. Along indicated route from nose towards the area near temple.



5. Along indicated route caring both sides of your nose wearing.

-5-