



## IPS Brightening Skin Corrector

User Manual

- Thank you for purchasing this beauty device.
- Please use under guidelines in this book and be well aware of the safety information listed. Keep this book for future reference.
- The usage of this device is premised on the basis of providing routine skin care. Any other ways of use and sale in connection with medical purpose are not supported.

### IPS Brightening Skin Corrector

Thank you for purchasing IPS Brightening Skin Corrector. Based on physiotherapy principles, this device is an exquisite and portable gadget which contains three most effective techniques which has been widely recognized in beauty industry: phototherapy, ultrasonic, and ion (galvanic). With it, you are able to make use of the harmless physical energy to promote the vitality of skin cells, to accelerate body circulation, and to strengthen the ability of oxygen supply and metabolism, so that the skin can regain young condition.

By using IPS Brightening Skin Corrector, we can easily achieve at home :

- Deeply cleansing;
- Hydrating and wrinkle smoothing;
- Tightening and lifting;
- Shrinking of puffiness and double-chin;
- Spots fading and skin tone brightening;
- Eye bags and dark circles removing.

### | CONTENTS |

1.Safety & Maintenance Precaution .....	02
2.Parts introduction .....	03
3.Photonic & Ultrasonic Anti-aging Device .....	04
4.Problems aimed at. ....	05
5.Working principles of the machine.....	05
6.Introduction of Ultrasonic .....	05
7.Introduction of Phototherapy .....	06
8.Introduction of Ionic care.....	06
9.Start using .....	07
10.Simple ways of testing allergic skin.....	07
11.Daily care .....	07
12.Special care of problematic skins .....	09
13.Spots fading, skin tone brightening .....	09
14.For sensitive skin .....	09
15.After use .....	09

### Product Details

Name: IPS Brightening Skin Corrector  
Input power: AC 100V ~ 240V  
Output power: DC 5V 1000 mA  
Dimensions (cm): 170\*46\*50mm  
Weight (g): 139  
Diameter of Ultrasound Head (mm): 28  
Ultrasonic frequency: 3MHz±5%


### Safety Instructions

- This device is designed to provide routine skin cares. Use for any other purpose is not supported.
- Please read this chapter carefully before use. Follow the guidelines while use and control the frequency of use under instruction of this manual. Please shorten the operating time for the first time of your use if your skin is very sensitive.
- Please stop using if you feel any uncomfortable. Consult your doctor right away.

### Do not use this device


- With any medical machines such as Cardiac Pacemaker, Oxygenator, electrocardiograph, etc.
- In any of circumstances below:
  - During pregnancy or menses
  - Have allergic dermatitis or hemorrhage
  - On Irritating skin or redness caused by sunburn
  - On skin ulcers or inflammations
  - On cosmetic operated body parts
  - Have telangiectasia or disordered liver function caused by long term uptake of steroid and hormone

### Operation & Maintenance Precautions

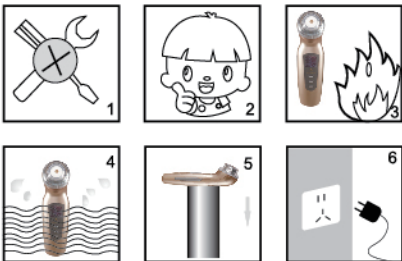

- Do not disassemble it by yourself
- Keep it away from children
- Keep it away from fire
- Keep it away from humidity and water
- Do not drop or collide it
- Unplug the electricity cord after charging

### Parts introduction

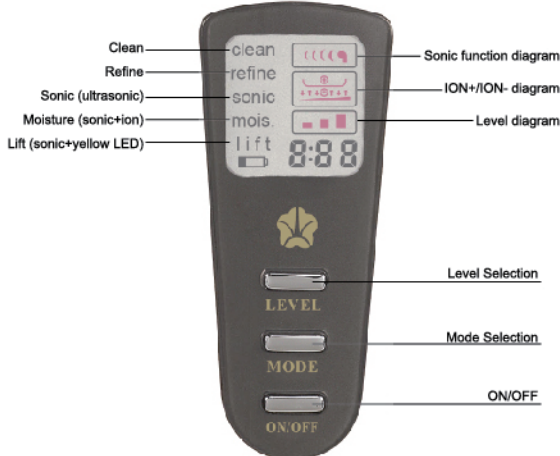

When a mode is selected, you can see from the display screen the following details:  
(1) The mode is underlined.  
(2) On the right side corresponding diagram shows.  
(3) Intensity level is defaulted to be low.  
(4) Suggested using time appears on the bottom of the screen.

- This Photonic & Ultrasonic Anti-aging Device includes the following parts:
1. Device (1)
  2. User manual (1)
  3. Charging set(1)
  4. Protective goggles (1 pair)

### Charging instruction

1. Before your first use, please charge the device for about 3 hours.
2. The charging type of this device is similar to that of a mobile phone. During the charging process, you can check the battery level from the battery icon on the screen. When it is fully charged, the battery icon will disappear.
3. When electricity is low, the screen will flash, so will the battery icon on the screen, accompanying with beeps. It turns off automatically after flashing for about 10 times.
4. When electricity falls low, charge it as soon as possible, otherwise the longtime low-voltage condition may cause damage to the battery.



### Skin concerns which IPS Brightening Skin Corrector is aiming at

This device is aiming at the following skin problems:

- Sag, looseness, wrinkles & lines
- Spots & blemishes
- Dull Skin tone
- Aging skin, slow metabolism, low nutrition-absorbing ability
- Dark circles & eye pouches

### Working principles of IPS Brightening Skin Corrector

#### Introduction of Ultrasonic function

Photonic & Ultrasonic Anti-aging Device has 3M Hz ultrasonic, which is the most effective and safest ultrasonic frequency recognized in the beauty industry. It can mildly penetrate into skin for cell massage and thus is particularly suitable for caring the skin around eyes. The high-speed micro-vibration brought by ultrasonic gives appropriate stimulation to skin cells, helping to improve the metabolism of cells, and activates skin vitality from inside.

### Introduction of phototherapy

#### Yellow light

The yellow photon under certain wavelength is especially suitable for caring sensitive skin. Yellow light can improve lymph system, strengthen function of muscle and immune system, and meanwhile it can increase the elasticity of capillaries, sooth irritation of skin, and brighten dull skin tone, fading spots and blemishes.

#### Introduction of Ionic care

Galvanic, ionic conducting is a technology that has been wildly accepted in beauty salons and hospitals. The difference in PH value between corneum cells and basal cells produces a weak electric field between the skin surface and the germinal layer, which blocks a lot of exotics including nutrition from going in deep skin. Galvanic, ionic function can also merge an electric field and when the two fields meet, the block is temporarily weakened and ionized nutrition is thus able to go into the germinal layer deeply. Working in the same function as a magnet with its positive and negative electrodes, galvanic and ionic treatment extracts the dirt and grease out of skin pores and penetrates nutrition deep into skin with an absorption rate many times higher than by hands or oral.

Moreover, the acid reaction caused by positive electrode can reduce the blood supply of skin thus relieve skin redness and irritations. It also can shrink the pores by strengthening the skin texture. This makes it greatly suitable for easily allergic skin and acne skin.

On the other side, the alkali reaction caused by negative electrode can stimulate the blood circulation and soften skin texture, thus enhance skin vitality and elasticity, and therefore greatly helps in aging skin and dryness.

### Start using

Before your first using, please do an allergic skin test. If allergy signs occur, you should stop using it immediately, and consult your doctor if necessary.

#### Allergic skin test

- Clean the wrist.
- Put your preferred skincare agent (gel/lotion) onto the cleansed area, turn on the device then choose to start the ultrasound function. Massage the area for 2-3 minutes, air-dry it and then observe for 48 hours. Stop using it immediately when you feel wrong with the skin.
- If no skin irritations or redness occur, it is ready to be used.



#### Before use

Lay the ultrasonic head flat with head upwards, drop a water onto the head, and press 'Sonic'. The water drop will vibrate when the ultrasound is on.

#### Caring Tips

- The surface of the probe should cling to the skin flatly.



- Apply gel to the skin and probe , increase the amount of gel when you feel the probe surface is not well-lubricated. (Suggested dosage of liquid is 4 to 5g at one time.)
- This device is designed to beep every 5 seconds as an indicator that it is working properly. It turns off automatically at the interval of 10 minutes. If you need to go on massaging, turn it on again.

#### Daily care

For daily facial care, the usage should not exceed three times per week, and each time of usage should not exceed 30 minutes.

#### Step 1 Deep cleansing

1. Remove your make-up and cleanse the face.
2. Pin the cotton pad onto the stainless steel probe.
3. Turn the device on, press MODE button to select CLEAN
4. Press LEVEL button to select a suitable intensity level for yourself.
5. Massage according to the picture below. (It takes around 5 seconds to massage from the chin to the root of ears, which also applies to other parts)
6. Cleanse the face and wipe the device clean.




#### Step 2 Lymph detox

- Apply gel onto your face as well as the probe.
- Turn the device on, press MODE button to select SONIC.
- Press LEVEL button to select a suitable intensity level, and do lymph detox according to the following picture.



#### Step 3 Nutrition conducting & lifting

- Apply skin care product onto your face and the probe.
- Turn the device on, press MODE to select MOISTURE or LIFT function.
- Press LEVEL to select a suitable intensity level.
- Massage according to the pictures below. (It takes 5 seconds to massage from the chin to the roots of ears, which also applies to other parts)

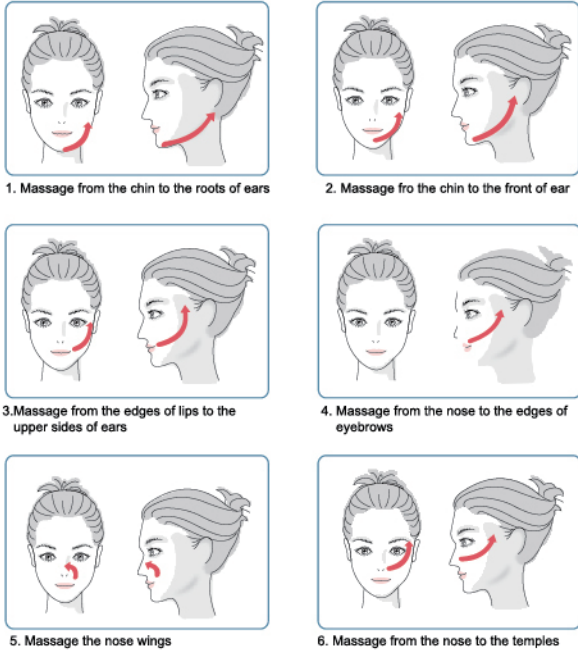

### Special care for problem skin

#### Step 1 Fading spots, correcting skin tone.

#### Step 2 Nutrition conducting

- Apply skin care product containing VC and/or other ingredients with spot-fading effect evenly onto the spots and the probe, or pin the cotton pad onto the probe directly.
- Turn the device on, press MODE button to select REFINE.
- Cling the probe to the spots' positions and massage, stay for 20 to 30 seconds respectively.

#### Step 3 Photonic care

Use yellow light (select LIFT mode and leave the skin untouched by the probe) to shine for 9 minutes. If you want more, reboot the LIFT mode.

#### Sensitive skin (use once every three days or per week)

- Select the single SONIC function, following the same steps as to those of daily care.
- Ion is not suggested to be used along.
- After caring, always remember to hydrate your skin well.

#### After use

- Turn off the device after use;
- Wipe the steel head gently with wet tissue.
- Do not clean it in water and/or with any chemical cleansing agent.
- Store the device clean and dry.